## Audífono inalámbrico binaural

Materiales: plástico ABS + silicona

Versión de Bluetooth: 5.0

Distancia de Bluetooth: 10 metros

Capacidad de la batería de auriculares: 35 mA

Capacidad de batería de caja de almacenamiento: 300 mA

Tiempo de juego: 4-5 horas


Cargos de cuadro de grabadora: 1-2 horas

Tiempo de espera: 120 horas

Nivel de NC: ☐ 25db Frecuencia: 2.4GHz

Tamaño de la caja: 73\*45\*25 mm



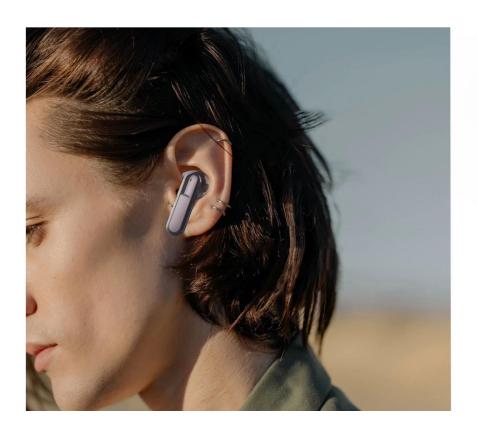






# Wireless Bluetooth transmission

Wireless Bluetooth 5.1, 10 meters distance transmission







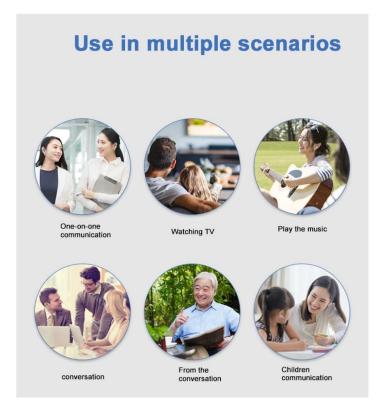







## **HOW TO USE**

Open the charging compartment, wait for the green/orange light to go out, and the headset prompts that the pairing is successful.

Open the charging compartment, wait for the green/orange light to go out, and the headset prompts that the pairing is successful. The earphones can be used individually or in pairs. In the state of playing music, touch the middle of the earphone handle with your hand, and the music sound of the music will change from small to large, and the volume can be controlled in a total of 16 levels

Pinch 1.5S to play music, then pinch 1.5S to pause the music

To answer the incoming call, tap 1S, and tap 1S to turn off the phone.

Hearing aid/assisted hearing mode, the 3S earphone has a prompt volume sound, indicating that the hearing aid auxiliary mode is turned on, and at the same time, the music can be played or paused without affecting the hearing aid function

In the hearing aid mode, the music volume can be adjusted.

You can use APP to pair earphones to adjust the effect of hearing aid magnification and adapt to different groups of people.

# TWS Binaural Wireless Hearing Aid




#### **IMPORTANT FUNCTIONS**


Wearing this hearing aid can significantly improve hearing

#### **OPERATION INTRODUCTION**

CHARGE WAREHOUSE

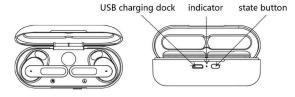

TWS EARPHONE

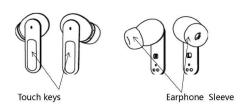

## **CHARGING COMPARTMENT FUNCTION**

NO	Function	operational approach	LED indication	
1	Charging the charging compartment	Connect lightning cable	When charging, the orange light flashes for a long time, and the orange light turns to green light for 3S after fully charged; Less than full power: the orange light keeps on for 8S;	
2	Stop charging	Press the trigger button to unplug the lightning cable	Green light is on last 3S	
3	Earphone in warehouse	Put the earphones in the compartment (after the earphones are in contact with the charging contacts)	Electricity> 10%, green light flash 1 time ; Electricity <=10%, orange light flash 1 time ;	
4	Headphone out of warehouse	Headphone release bin (after the headphone leaves the charging contacts)	Electricity> 10%, green lightflash 1 time ; Electricity <=10%, orange light flash 1 time ;	
5	low battery	When the power level of the trigger button is 0, remind	Orange light last flish, 3s;	
6	Charging compartment battery check	Short press button	Electricity> 25%, green light flash 1 time ; Electricity <=25%, orange light flash 1 time ;	
7	Uncap	Open the charging compartment cover	Electricity> 10%, the green light turns on for 3 seconds and then goes out; Electricity <=10%, range light Flash 3 times;	
8	Close the lid	Close the charging compartment cover	Electricity> 10%, the green light turns on for 3 seconds and then goes out; Electricity <=10%, the orange light turns on for 3s and goes out;	
9	In-house matching	Open the cover, the green light/orange light is on and the headset is paired in the compartment	Pairing successfully: the headset will prompt the pairing success message ONE TWO; Pairing failed: no prompt;	
10	Full battery prompt	After fully charged	When it reaches full power: the green light is always on; When the lightning line is connected: When the bastery is not full, press the button shortly, the orange light flashes for a long time; Fully charged, short press the button, the green light flash 1 time.	

#### PROMPT TONE

Key pattern	Function	functional operation	LED indicate	voice prompt
	Power on	open a position	The green light goes out	POWER ON
	Power off	Warehousing	Green light out	POWER OFF
	Match	Warehouse Entry/Exit	Green light out	ONE TWO
	Hearing Aid Mode	Default	NA	
MFB	Pure Bluetooth mode	Touch &Hold 3s	NA	The beep goes on twice
	Answerin g calls	Short Press	NA	One beep
	hang up	Long press 1s/ double click	NA	One beep
	Decline calls	Long press 1 NA	NA	One beep
	play music	Long press 1.5s	NA	One beep
	Pause	Long press 1.5s	NA	One beep
V+/V-	volume plus	Double click for 0.5s / hearing aid volume plus one level, a total of 6 levels, the volume is cycled from high to low	NA	Beep

## **HEARING AID SPECIFICATIONS**

	sensitivity of a pickup	-38dB	
	Pickup distance	3-5m	
Hearing aid specificatio	Pickup noise suppression	-60db	
ns	Howling suppression ratio	-40db	
	Pickup gain	-30db	
	Pickup frequency	125HZ-10KHZ	

## **BLUETOOTH SPECLFICATIONS**

	DSP	QC9864BT	
	Bluetooth version	5.1	
Bluetooth	Fraquency range	2.402GHz-2.480GHz	
specifications	ransmission distance (unobstructed)	≤15meter	
	Track system	Stereo	
	Speaker Impedance	32Ω/Φ10mm	
	Audio Power Output(Max)	135mW	
Output	Headphones Sensitivity	120dB	
specifications	Fraquency range	20HZ-20KHZ	
	distortion factor	<1%	
	SNR	≥90dB	
	Depth of call noise reduction	adjustable0-30dB	

## PREGUNTAS MÁS FRECUENTES:

## Acerca de la compañía:

- ¿Cuánto tiempo hace que su empresa trata con negocios electrónicos? ¿Estás fabricando? Jaskey Limited es una empresa de diseño, desarrollo de desarrollo durante 12 años y nos especializamos en productos electrónicos de moda modernos con valor agregado al mercado electrónico de consumo.
- ¿Qué auditoría de fábrica hiciste?

Nuestra fabricación ya realizó una auditoría de fábrica de Sedex 4P, BSCI, ISO9001 con un sistema de control de calidad profesional.

- ¿Qué mercado es su principal mercado comercial?
   Nuestro principal mercado comercial es Europa, EE. UU., Asia, América del Sur, etc.
- ¿Su empresa tiene su propia marca o no?
   Tenemos nuestra propia marca Sawolol para nuestra sección de audio inteligente de diseño patentado.
- ¿Qué plazo de pago para su empresa?

Normalmente es T/T por adelantado (depósito del 30% antes de la confirmación del pedido y el saldo antes del envío) para el pedido de masa o LC a la vista para un pedido de gran cantidad; Para orden de muestra con una pequeña cantidad, está bien transferir el pago a través de PayPal, Western Union o Alipay también.

- ¿A qué feria asistirá su empresa cada año?

Normalmente asistiremos a HKTDC Hong Kong Electronics Fair en la edición de primavera y otoño y la feria de CES, etc., mantenga el contacto con nosotros para obtener nuestra información actualizada a través de nuestro sitio web o nuestro correo info@jaskeyworld.com , iGracias!

# Sobre productos:

- ¿Puedes hacer productos OEM/ODM?

Tenemos una gran experiencia en productos de audio personalizados OEM/ODM con equipo de I + D, equipo de diseño de productos que puede proporcionar 3 nuevos productos de diseño trimestralmente con diseño único y función de alta calidad a nuestro socio comercial.

- ¿Qué certificado para su desinfectante automático? Nuestros productos se cumplen con los productos electrónicos que prueban el estándar de CE, Rosh, FCC, BQB, etc.
- ¿Cuál es el tiempo de entrega para el orden normal?
   Tiempo de muestra: 3-7 trabajando para muestras normales.
   Tiempo de pedido de masa: 25-35 días después de la aprobación de la muestra para pedidos diarios.
- ¿Qué paquete para el desinfectante automático? Normalmente es nuestra caja de regalo estándar y está bien hacer un paquete personalizado según los requisitos.
- ¿Cuánto tiempo para la garantía del producto?

La garantía de los productos es un año a partir de la fecha de envío de los productos bajo uso normal, excluyendo baterías y otros accesorios prescindibles.

Por favor, comuníquese con nosotros libremente para obtener cualquier consulta a través de nuestro correo de servicio al cliente info@jaskeyworld.com, igracias! iLo apoyaremos para nuestra próxima cooperación comercial, igracias!

Este es nuestros productos de venta en caliente:

Aurel inalámbrico de cancelación de ruido activo AEP-0219

Auriculares TWS de cancelación de ruido activo Bluetooth v5.1 Sport inalámbrico y auriculares Gaming AEP-0244

Auriculares de conducción de huesos de Openear HEP-0170